

# anti- AMPK beta 1 antibody

### **Product Information**

Catalog No.:	FNab00377
Size:	100µg
Form:	liquid
Purification:	Immunogen affinity purified
Purity:	$\geq$ 95% as determined by SDS-PAGE
Host:	Rabbit
Clonality:	polyclonal
Clone ID:	None
IsoType:	IgG
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3, -20°C for 12 months (Avoid repeated freeze / thaw cycles.)

### Background

The protein encoded by this gene is a regulatory subunit of the AMP-activated protein kinase (AMPK). AMPK is a heterotrimer consisting of an alpha catalytic subunit, and non-catalytic beta and gamma subunits. AMPK is an important energy-sensing enzyme that monitors cellular energy status. In response to cellular metabolic stresses, AMPK is activated, and thus phosphorylates and inactivates acetyl-CoA carboxylase (ACC) and beta-hydroxy beta-methylglutaryl-CoA reductase (HMGCR), key enzymes involved in regulating de novo biosynthesis of fatty acid and cholesterol. This subunit may be a positive regulator of AMPK activity. The myristoylation and phosphorylation of this subunit have been shown to affect the enzyme activity and cellular localization of AMPK. This subunit may also serve as an adaptor molecule mediating the association of the AMPK complex.

#### **Immunogen information**

Immunogen:	protein kinase, AMP-activated, beta 1 non-catalytic subunit
Synonyms:	AMPK, AMPK beta 1, AMPK subunit beta 1, AMPKb, HAMPKb, PRKAB1
Observed MW:	37 kDa
Uniprot ID :	Q9Y478

#### Wuhan Fine Biotech Co., Ltd.


B9 Bld, High-Tech Medical Devices Park, No. 818 Gaoxin Ave.East Lake High-Tech Development Zone.Wuhan, Hubei, China(430206)

1




## Application

Reactivity:Human, Mouse, RatTested Application:ELISA, WB, IHCRecommended dilution:WB: 1:500 - 1:2000; IHC: 1:50 - 1:200Image:

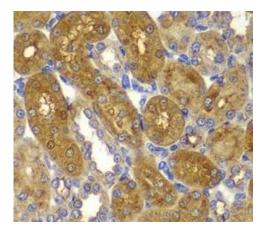

Immunohistochemistry of paraffin-embedded rat kidney using FNab00377(PRKAB1 antibody) at dilution of 1:100

mouse liver tissue were subjected to SDS PAGE followed by western blot with FNab00377( PRKAB1 Antibody) at dilution of 1:1000

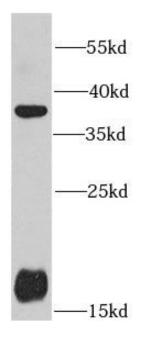

Wuhan Fine Biotech Co., Ltd.

B9 Bld, High-Tech Medical Devices Park, No. 818 Gaoxin Ave.East Lake High-Tech Development Zone.Wuhan, Hubei, China(430206)

Tel :( 0086)027-87384275 Fax: (0086)027-87800889

www.fn-test.com