

BAIAP3 antibody

Product Information

Catalog No.:	FNab00795
Size:	100µg
Form:	liquid
Purification:	Immunogen affinity purified
Purity:	≥95% as determined by SDS-PAGE
Host:	Rabbit
Clonality:	polyclonal
Clone ID:	None
IsoType:	IgG
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3, -20°C for 12 months(Avoid repeated freeze / thaw cycles.)

Background


None

Immunogen information


Immunogen:	BAI1-associated protein 3
Synonyms:	BAI1-associated protein 3 (BAP3) Brain-specific angiogenesis inhibitor I-associated protein 3 BAIAP3 KIAA0734
Observed MW:	130 kDa
Uniprot ID :	O94812

Application

Reactivity:	Human, Mouse, Rat
Tested Application:	ELISA, WB, IF
Recommended dilution:	WB: 1:200-1:1000; IF: 1:10-1:100
Image:	

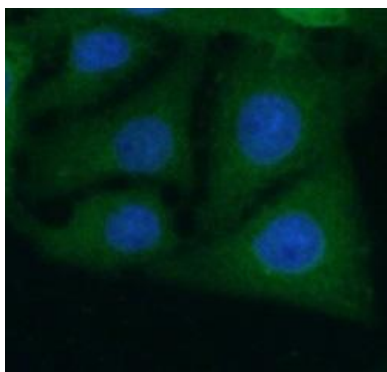

Immunofluorescent analysis of SH-SY5Y cells using FNab00795(BAIAP3 antibody) at dilution of 1:25 and Alexa Fluor 488-conjugated Goat Anti-Rabbit IgG(H+L)

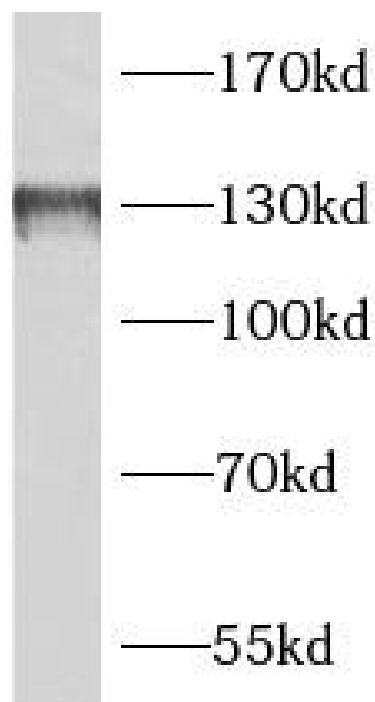

mouse brain tissue were subjected to SDS PAGE followed by western blot with FNab00795(BAIAP3 Antibody) at dilution of 1:300