

SRC antibody

Product Information


Catalog No.:	FNab02028
Size:	100µg
Form:	liquid
Purification:	Protein A+G purification
Purity:	≥95% as determined by SDS-PAGE
Host:	Mouse
Clonality:	monoclonal
Clone ID:	9H2
IsoType:	IgG1
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3, -20°C for 12 months(Avoid repeated freeze / thaw cycles.)

Immunogen information


Immunogen:	v-src sarcoma(Schmidt-Ruppin A-2) viral oncogene homolog(avian)
Synonyms:	Proto-oncogene tyrosine-protein kinase Src Proto-oncogene c-Src pp60c-src (p60-Src) SRC SRC1
Observed MW:	60 kDa
Uniprot ID :	P12931

Application

Reactivity:	Human, Hamster
Tested Application:	ELISA, WB, IHC, IF
Recommended dilution:	WB: 1:1000-1:8000; IHC: 1:100-1:400; IF: 1:50-1:500
Image:	

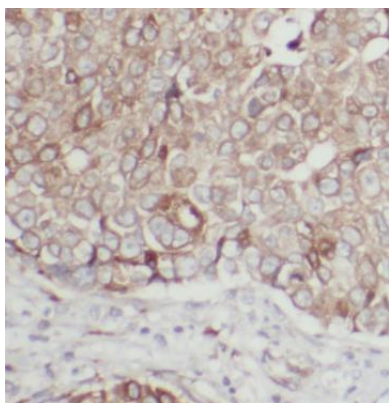

Immunohistochemistry of paraffin-embedded human breast cancer tissue slide using FNab02028(c-SRC Antibody) at dilution of 1:200 heat mediated antigen retrieved with Tris-EDTA buffer(pH9).

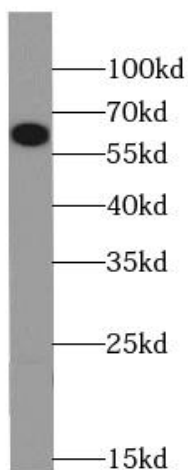

Western blot of c-SRC in various cell lines with FNab02028 at dilution of 1:4000. FNab02028