

## GABARAPL1 antibody

### Product Information

Catalog No.:	FNab03275
Size:	100µg
Form:	liquid
Purification:	Immunogen affinity purified
Purity:	≥95% as determined by SDS-PAGE
Host:	Rabbit
Clonality:	polyclonal
Clone ID:	None
IsoType:	IgG
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3, -20°C for 12 months(Avoid repeated freeze / thaw cycles.)

### Background

Ubiquitin-like modifier that increases cell-surface expression of kappa-type opioid receptor through facilitating anterograde intracellular trafficking of the receptor. Involved in formation of autophagosomal vacuoles. Whereas LC3s are involved in elongation of the phagophore membrane, the GABARAP/GATE-16 subfamily is essential for a later stage in autophagosome maturation.


### Immunogen information

Immunogen:	GABA(A) receptor-associated protein like 1
Synonyms:	Gamma-aminobutyric acid receptor-associated protein-like 1 Early estrogen-regulated protein GABA(A) receptor-associated protein-like 1 Glandular epithelial cell protein 1 (GEC-1) GABARAPL1 GEC1
Observed MW:	17 kDa
Uniprot ID :	Q9H0R8


### Application

Reactivity:	Human, Mouse, Rat
Tested Application:	ELISA, WB, IHC
Recommended dilution:	WB: 1:500-1:2000; IHC: 1:20-1:200

Image:

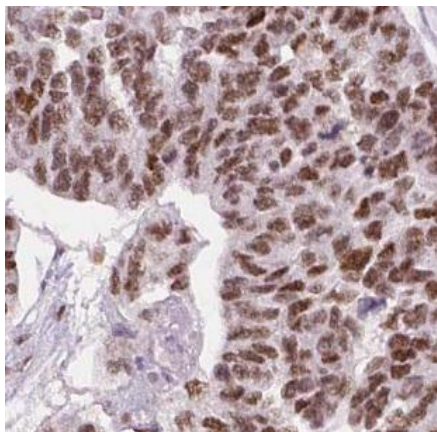

Immunohistochemistry of paraffin-embedded human ovary tumor using FNab03275(GABARAPL1 antibody) at dilution of 1:100

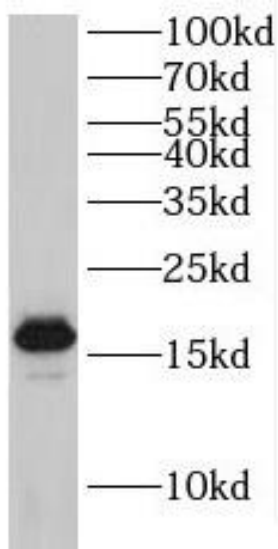

mouse kidney tissue were subjected to SDS PAGE followed by western blot with FNab03275(GABARAPL1 antibody) at dilution of 1:500