

GPX4 antibody

Product Information

Catalog No.:	FNab03622
Size:	100µg
Form:	liquid
Purification:	Immunogen affinity purified
Purity:	≥95% as determined by SDS-PAGE
Host:	Rabbit
Clonality:	polyclonal
Clone ID:	None
IsoType:	IgG
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3, -20°C for 12 months (Avoid repeated freeze / thaw cycles.)

Background

The protein encoded by this gene belongs to the glutathione peroxidase family, members of which catalyze the reduction of hydrogen peroxide, organic hydroperoxides and lipid hydroperoxides, and thereby protect cells against oxidative damage. Several isozymes of this gene family exist in vertebrates, which vary in cellular location and substrate specificity. This isozyme has a high preference for lipid hydroperoxides and protects cells against membrane lipid peroxidation and cell death. It is also required for normal sperm development; thus, it has been identified as a 'moonlighting' protein because of its ability to serve dual functions as a peroxidase, as well as a structural protein in mature spermatozoa. Mutations in this gene are associated with Sedaghatian type of spondylometaphyseal dysplasia (SMDS). This isozyme is also a selenoprotein, containing the rare amino acid selenocysteine (Sec) at its active site. Sec is encoded by the UGA codon, which normally signals translation termination. The 3' UTRs of selenoprotein mRNAs contain a conserved stem-loop structure, designated the Sec insertion sequence (SECIS) element, that is necessary for the recognition of UGA as a Sec codon, rather than as a stop signal. Alternatively spliced transcript variants have been found for this gene.

Immunogen information

Immunogen:	glutathione peroxidase 4 (phospholipid hydroperoxidase)
Synonyms:	Phospholipid hydroperoxide glutathione peroxidase GPX4 (PHGPx) Glutathione peroxidase 4 (GPx-4, GSHPx-4) GPX4
Observed MW:	25 kDa
Uniprot ID :	P36969


Application

Reactivity: Human, Mouse, Rat


Tested Application: ELISA, WB, IHC, IF

Recommended dilution: WB: 1:500 - 1:2000; IHC: 1:50 - 1:200; IF: 1:50 - 1:200

Image:

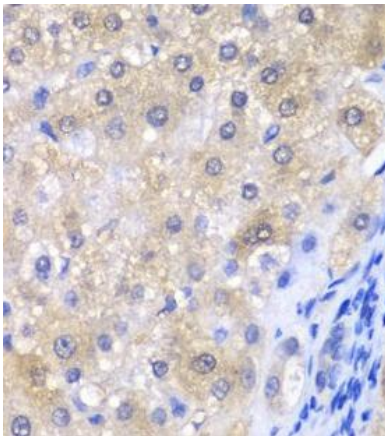

Immunohistochemistry of paraffin-embedded human liver injury using FNab03622(GPX4 antibody) at dilution of 1:50

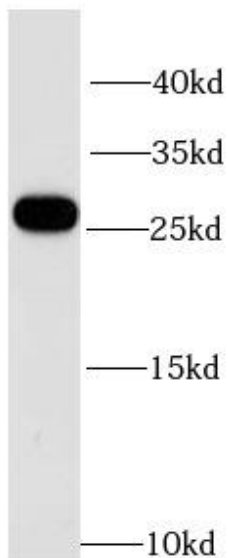

mouse kidney tissue were subjected to SDS PAGE followed by western blot with FNab03622(GPX4 antibody) at dilution of 1:1000