

TRIM28 antibody

Product Information

Catalog No.:	FNab04466
Size:	100µg
Form:	liquid
Purification:	Immunogen affinity purified
Purity:	≥95% as determined by SDS-PAGE
Host:	Rabbit
Clonality:	polyclonal
Clone ID:	None
IsoType:	IgG
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3, -20°C for 12 months(Avoid repeated freeze / thaw cycles.)

Background


None

Immunogen information


Immunogen:	tripartite motif-containing 28
Synonyms:	Transcription intermediary factor 1-beta (TIF1-beta) E3 SUMO-protein ligase TRIM28 KRAB-associated protein 1 (KAP-1) KRAB-interacting protein 1 (KRIP-1) Nuclear corepressor KAP-1 RING finger protein 96 RING-type E3 ubiquitin transferase TIF1-beta Tripartite motif-containing protein 28 TRIM28 KAP1 RNF96 TIF1B
Observed MW:	100 kDa
Uniprot ID :	Q13263

Application


Reactivity:	Human, Mouse, Rat
Tested Application:	ELISA, WB, IHC, IF, FC, IP
Recommended dilution:	WB: 1:500-1:2000; IP: 1:200-1:1000; IHC: 1:20-1:200; IF: 1:20-1:200
Image:	

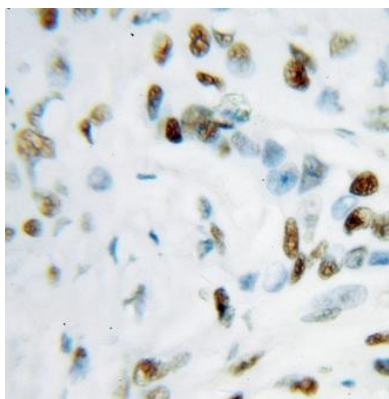

Immunohistochemistry of paraffin-embedded human colon cancer using FNab04466(KAP1 antibody) at dilution of 1:100

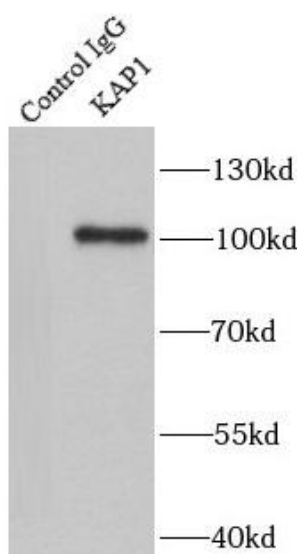

IP Result of anti-KAP1 (IP:FNab04466, 4ug; Detection:FNab04466 1:1000) with HeLa cells lysate 1200ug.

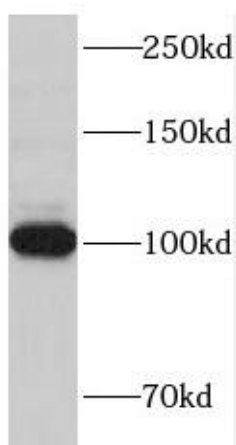

HeLa cells were subjected to SDS PAGE followed by western blot with FNab04466(KAP1 antibody) at dilution of 1:3000