

## LY6D antibody

### Product Information

Catalog No.:	FNab04894
Size:	100µg
Form:	liquid
Purification:	Immunogen affinity purified
Purity:	≥95% as determined by SDS-PAGE
Host:	Rabbit
Clonality:	polyclonal
Clone ID:	None
IsoType:	IgG
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3, -20°C for 12 months(Avoid repeated freeze / thaw cycles.)

### Background


May act as a specification marker at earliest stage specification of lymphocytes between B-and T-cell development. Marks the earliest stage of B-cell specification.

### Immunogen information


Immunogen:	lymphocyte antigen 6 complex, locus D
Synonyms:	Lymphocyte antigen 6D (Ly-6D) E48 antigen LY6D E48
Observed MW:	13-14 kDa
Uniprot ID :	Q14210

### Application

Reactivity:	Human, Mouse
Tested Application:	ELISA, WB, IHC
Recommended dilution:	WB: 1:200-1:1000; IHC: 1:20-1:200
Image:	

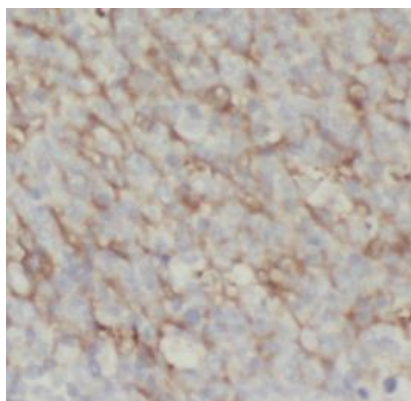

Immunohistochemistry of paraffin-embedded human tonsillitis tissue slide using FNab04894(LY6D Antibody) at dilution of 1:50

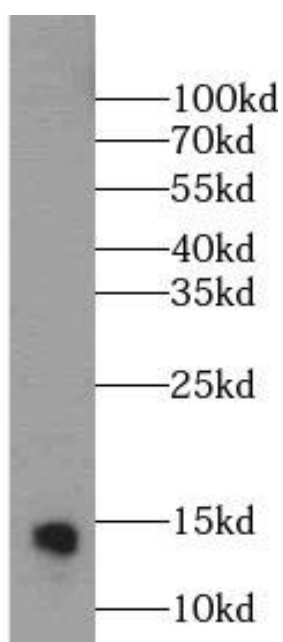

mouse spleen tissue were subjected to SDS PAGE followed by western blot with FNab04894(LY6D Antibody) at dilution of 1:300