

NUDT9 antibody

Product Information

Catalog No.:	FNab05908
Size:	100µg
Form:	liquid
Purification:	Immunogen affinity purified
Purity:	≥95% as determined by SDS-PAGE
Host:	Rabbit
Clonality:	polyclonal
Clone ID:	None
IsoType:	IgG
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3, -20°C for 12 months(Avoid repeated freeze / thaw cycles.)

Background


Hydrolyzes ADP-ribose(ADPR) to AMP and ribose 5'-phosphate.

Immunogen information


Immunogen:	nudix(nucleoside diphosphate linked moiety X)-type motif 9
Synonyms:	ADP-ribose pyrophosphatase, mitochondrial ADP-ribose diphosphatase ADP-ribose phosphohydrolase Adenosine diphosphoribose pyrophosphatase (ADPR-PPase) Nucleoside diphosphate-linked moiety X motif 9 (Nudix motif 9) NUDT9 NUDT10
Observed MW:	39 kDa
Uniprot ID :	Q9BW91

Application

Reactivity:	Human, Mouse, Rat
Tested Application:	ELISA, WB, IHC
Recommended dilution:	WB: 1:500-1:2000; IHC: 1:20-1:200
Image:	

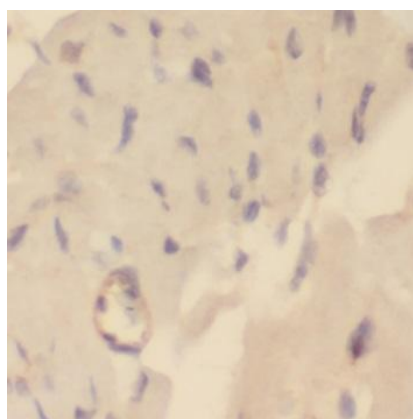

Immunohistochemistry of paraffin-embedded human heart using FNab05908(NUDT9 antibody) at dilution of 1:50

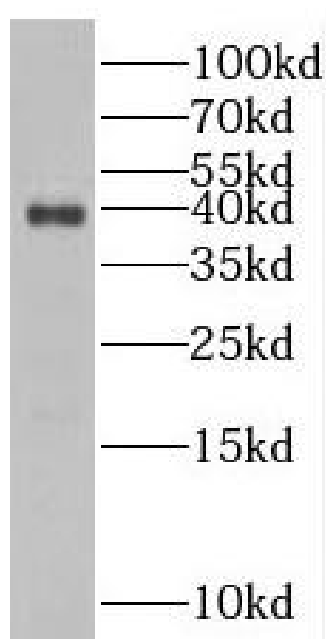

mouse brain tissue were subjected to SDS PAGE followed by western blot with FNab05908(NUDT9 antibody) at dilution of 1:1000