

## PITPNA antibody

### Product Information

Catalog No.:	FNab06470
Size:	100µg
Form:	liquid
Purification:	Immunogen affinity purified
Purity:	≥95% as determined by SDS-PAGE
Host:	Rabbit
Clonality:	polyclonal
Clone ID:	None
IsoType:	IgG
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3, -20°C for 12 months(Avoid repeated freeze / thaw cycles.)

### Background


Catalyzes the transfer of PtdIns and phosphatidylcholine between membranes.

### Immunogen information


Immunogen:	phosphatidylinositol transfer protein, alpha
Synonyms:	Phosphatidylinositol transfer protein alpha isoform (PI-TP-alpha, PtdIns transfer protein alpha, PtdInsTP alpha) PITPNA PITPN
Observed MW:	32 kDa
Uniprot ID :	Q00169

### Application

Reactivity:	Human, Mouse, Rat
Tested Application:	ELISA, WB, IHC
Recommended dilution:	WB: 1:500-1:2000; IHC: 1:20-1:200
Image:	

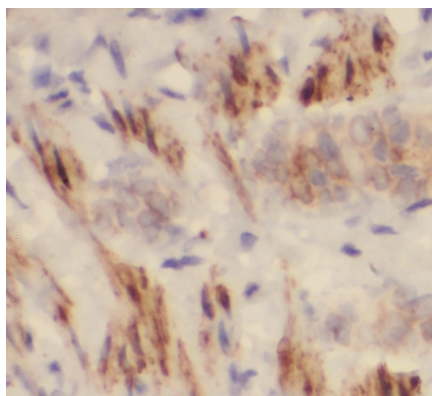

Immunohistochemistry of paraffin-embedded human lung using FNab06470(PITPNA antibody) at dilution of 1:100

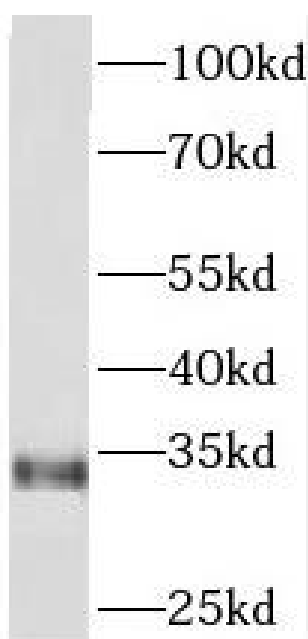

HEK-293 cells were subjected to SDS PAGE followed by western blot with FNab06470(PITPNA antibody) at dilution of 1:800