

## PYGO2 antibody

### Product Information

Catalog No.:	FNab06978
Size:	100µg
Form:	liquid
Purification:	Immunogen affinity purified
Purity:	≥95% as determined by SDS-PAGE
Host:	Rabbit
Clonality:	polyclonal
Clone ID:	None
IsoType:	IgG
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3, -20°C for 12 months(Avoid repeated freeze / thaw cycles.)

### Background


Involved in signal transduction through the Wnt pathway.

### Immunogen information


Immunogen:	pygopus homolog 2(Drosophila)
Synonyms:	Pygopus homolog 2 PYGO2
Observed MW:	50 kDa
Uniprot ID :	Q9BRQ0

### Application

Reactivity:	Human, Mouse, Rat
Tested Application:	ELISA, IHC, WB
Recommended dilution:	WB: 1:500-1:2000; IHC: 1:20-1:200
Image:	

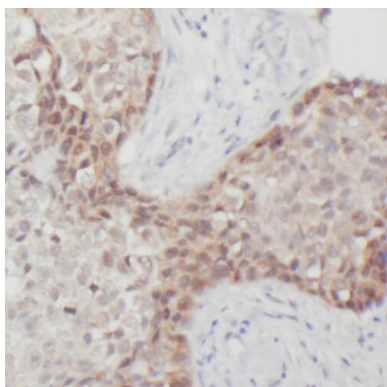

Immunohistochemistry of paraffin-embedded human ovary tumor tissue slide using FNab06978(PYGO2 Antibody) at dilution of 1:200 heat mediated antigen retrieved with Tris-EDTA buffer(pH9).

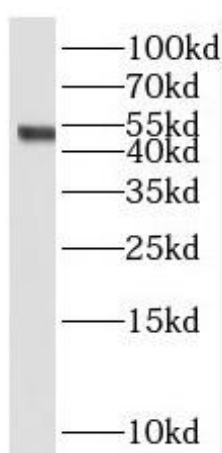

MDA-MB-453s cells were subjected to SDS PAGE followed by western blot with FNab06978(PYGO2 Antibody) at dilution of 1:600