

RBM24 antibody

Product Information

Catalog No.:	FNab07170
Size:	100µg
Form:	liquid
Purification:	Immunogen affinity purified
Purity:	≥95% as determined by SDS-PAGE
Host:	Rabbit
Clonality:	polyclonal
Clone ID:	None
IsoType:	IgG
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3, -20°C for 12 months(Avoid repeated freeze / thaw cycles.)

Background


Plays a role in myogenic differentiation by regulating MYOG levels. Binds to the 3'-UTR of MYOG mRNA and regulates its stability.

Immunogen information

Immunogen:	RNA binding motif protein 24
Synonyms:	RNA-binding protein 24 RNA-binding motif protein 24 RNA-binding region-containing protein 6 RBM24 RNPC6
Observed MW:	None
Uniprot ID :	Q9BX46

Application

Reactivity:	Human, Mouse, Rat
Tested Application:	ELISA, IHC
Recommended dilution:	IHC: 1:20-1:200
Image:	

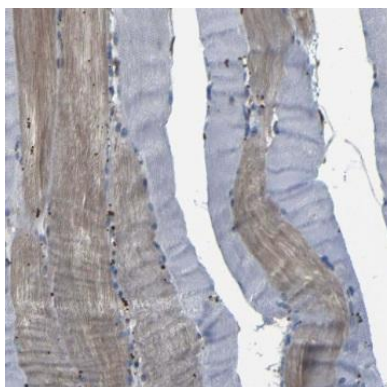

Immunohistochemistry of paraffin-embedded human skeletal muscle tissue using FNab07170(RBM24 antibody) at dilution of 1:50