

RSU1 antibody

Product Information

Catalog No.:	FNab07514
Size:	100µg
Form:	liquid
Purification:	Immunogen affinity purified
Purity:	≥95% as determined by SDS-PAGE
Host:	Rabbit
Clonality:	polyclonal
Clone ID:	None
IsoType:	IgG
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3, -20°C for 12 months(Avoid repeated freeze / thaw cycles.)

Background


Potentially plays a role in the Ras signal transduction pathway. Capable of suppressing v-Ras transformation in vitro.

Immunogen information


Immunogen:	Ras suppressor protein 1
Synonyms:	Ras suppressor protein 1 (RSP-1, Rsu-1) RSU1 RSP1
Observed MW:	32 kDa
Uniprot ID :	Q15404

Application

Reactivity:	Human, Mouse
Tested Application:	ELISA, WB, IP
Recommended dilution:	WB: 1:500-1:2000; IP: 1:200-1:2000
Image:	

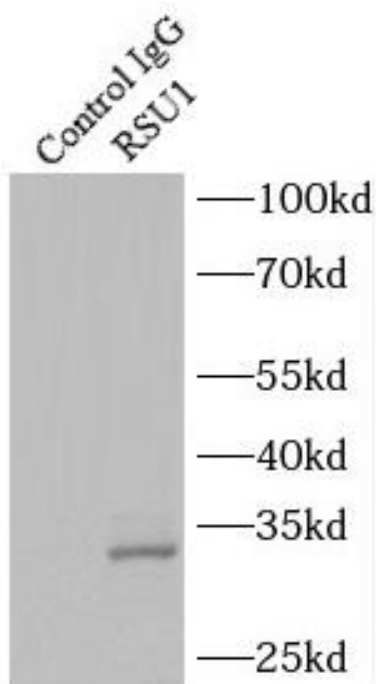

IP Result of anti-RSU1 (IP:FNab07514, 4ug;
Detection:FNab07514 1:300) with human placenta
tissue lysate 1520ug.

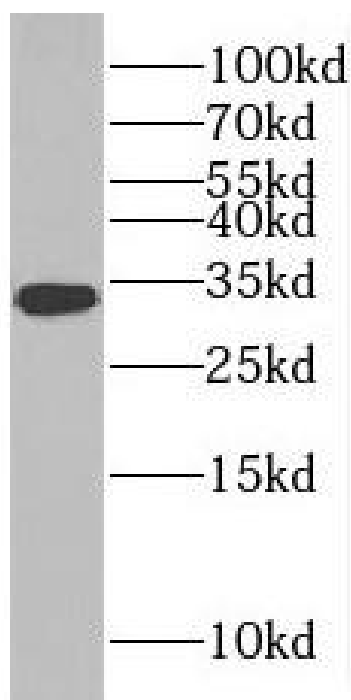

mouse placenta tissue were subjected to SDS
PAGE followed by western blot with
FNab07514(RSU1 antibody) at dilution of 1:300