

## STK25 antibody

### Product Information

Catalog No.:	FNab08328
Size:	100µg
Form:	liquid
Purification:	Immunogen affinity purified
Purity:	≥95% as determined by SDS-PAGE
Host:	Rabbit
Clonality:	polyclonal
Clone ID:	None
IsoType:	IgG
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3, -20°C for 12 months(Avoid repeated freeze / thaw cycles.)

### Background


Oxidant stress-activated serine/threonine kinase that may play a role in the response to environmental stress. Targets to the Golgi apparatus where it appears to regulate protein transport events, cell adhesion, and polarity complexes important for cell migration.

### Immunogen information

Immunogen:	serine/threonine kinase 25(STE20 homolog, yeast)
Synonyms:	Serine/threonine-protein kinase 25 Ste20-like kinase Sterile 20/oxidant stress-response kinase 1 (SOK-1, Ste20/oxidant stress response kinase 1) STK25 SOK1 YSK1
Observed MW:	48 kDa
Uniprot ID :	O00506

### Application

Reactivity:	Human, Mouse
Tested Application:	ELISA, WB
Recommended dilution:	WB: 1:500-1:2000
Image:	

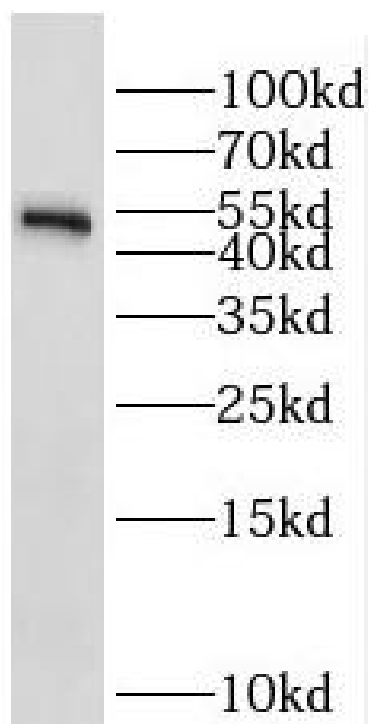

Neuro-2a cells were subjected to SDS PAGE followed by western blot with FNab08328(STK25 antibody) at dilution of 1:800