

anti- TBX2 antibody

Product Information

Catalog No.:	FNab08536
Size:	100µg
Form:	liquid
Purification:	Immunogen affinity purified
Purity:	≥95% as determined by SDS-PAGE
Host:	Rabbit
Clonality:	polyclonal
Clone ID:	None
IsoType:	IgG
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3, -20°C for 12 months (Avoid repeated freeze / thaw cycles.)

Background

Involved in the transcriptional regulation of genes required for mesoderm differentiation. Probably plays a role in limb pattern formation. Acts as a negative regulator of PML function in cellular senescence. May be required for cardiac atrioventricular canal formation.

Immunogen information

Immunogen:	T-box 2
Synonyms:	TBX2, T box 2, T box protein 2, FLJ10169
Observed MW:	74 kDa
UniprotID :	Q13207

Application


Reactivity:	Human, Mouse, Rat
Tested Application:	ELISA, WB, IHC
Recommended dilution:	WB: 1:200-1:2000; IHC: 1:20-1:200
Image:	

Wuhan Fine Biotech Co., Ltd.


B9 Bld, High-Tech Medical Devices Park, No. 818 Gaoxin Ave. East Lake High-Tech Development Zone. Wuhan, Hubei, China (430206)

Tel : (0086)027-87384275

Fax: (0086)027-87800889 www.fn-test.com

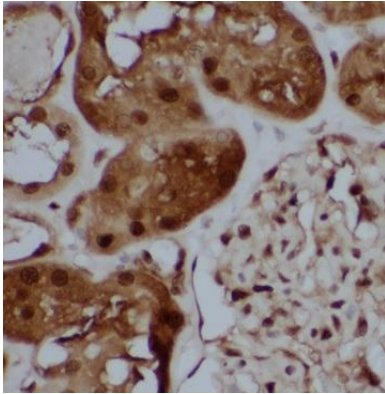

Immunohistochemistry of paraffin-embedded human kidney using FNab08536(TBX2 antibody) at dilution of 1:100

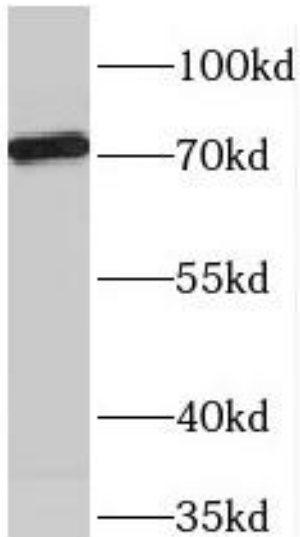

mouse lung tissue were subjected to SDS PAGE followed by western blot with FNab08536(TBX2 Antibody) at dilution of 1:300