

TMEM9B antibody

Product Information

Catalog No.:	FNab08794
Size:	100µg
Form:	liquid
Purification:	Immunogen affinity purified
Purity:	≥95% as determined by SDS-PAGE
Host:	Rabbit
Clonality:	polyclonal
Clone ID:	None
IsoType:	IgG
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3, -20°C for 12 months(Avoid repeated freeze / thaw cycles.)

Background


Enhances production of proinflammatory cytokines induced by TNF, IL1B, and TLR ligands.
Has a role in TNF activation of both the NF-kappaB and MAPK pathways.

Immunogen information


Immunogen:	TMEM9 domain family, member B
Synonyms:	Transmembrane protein 9B TMEM9B C11orf15
Observed MW:	23, 35 kDa
Uniprot ID :	Q9NQ34

Application

Reactivity:	Human, Rat, Mouse
Tested Application:	ELISA, WB, IHC
Recommended dilution:	WB: 1:500-1:2000; IHC: 1:20-1:200
Image:	

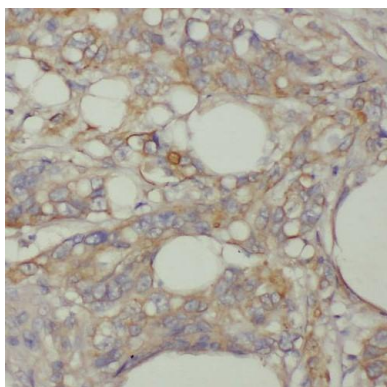

Immunohistochemistry of paraffin-embedded human colon cancer slide using FNab08794 (TMEM9B Antibody) at dilution of 1:50

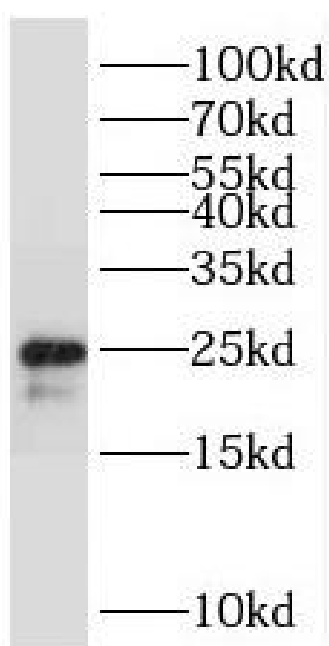

rat liver tissue were subjected to SDS PAGE followed by western blot with FNab08794 (TMEM9B Antibody) at dilution of 1:1000