

UBAC1 antibody

Product Information

Catalog No.:	FNab09147
Size:	100µg
Form:	liquid
Purification:	Immunogen affinity purified
Purity:	≥95% as determined by SDS-PAGE
Host:	Rabbit
Clonality:	polyclonal
Clone ID:	None
IsoType:	IgG
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3, -20°C for 12 months(Avoid repeated freeze / thaw cycles.)

Background


Non-catalytic subunit of the KPC complex that acts as E3 ubiquitin-protein ligase. Required for poly-ubiquitination and proteasome-mediated degradation of CDKN1B during G1 phase of the cell cycle.

Immunogen information


Immunogen:	UBA domain containing 1
Synonyms:	Ubiquitin-associated domain-containing protein 1 (UBA domain-containing protein 1) Glioblastoma cell differentiation-related protein 1 Kip1 ubiquitination-promoting complex protein 2 UBAC1 GBDR1 KPC2 UBADC1
Observed MW:	45 kDa
Uniprot ID :	Q9BSL1

Application

Reactivity:	Human, Mouse, Rat
Tested Application:	ELISA, WB, IHC
Recommended dilution:	WB: 1:500-1:2000; IHC: 1:20-1:200
Image:	

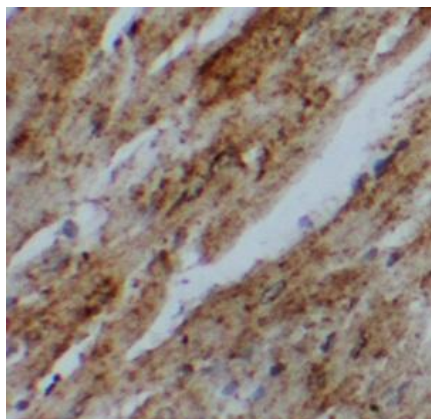

Immunohistochemistry of paraffin-embedded human heart using FNab09147(UBAC1 antibody) at dilution of 1:100

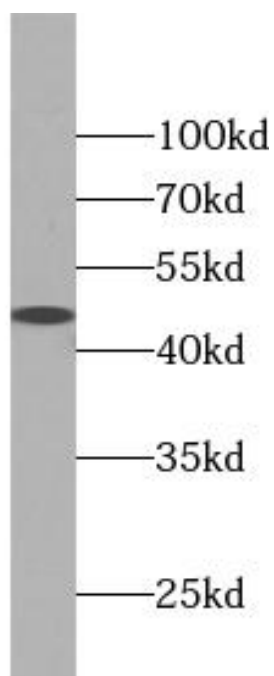

human brain tissue were subjected to SDS PAGE followed by western blot with FNab09147(UBAC1 antibody) at dilution of 1:500