

VBP1 antibody

Product Information

Catalog No.:	FNab09379
Size:	100µg
Form:	liquid
Purification:	Immunogen affinity purified
Purity:	≥95% as determined by SDS-PAGE
Host:	Rabbit
Clonality:	polyclonal
Clone ID:	None
IsoType:	IgG
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3, -20°C for 12 months(Avoid repeated freeze / thaw cycles.)

Background


Binds specifically to cytosolic chaperonin(c-CPN) and transfers target proteins to it. Binds to nascent polypeptide chain and promotes folding in an environment in which there are many competing pathways for nonnative proteins.

Immunogen information


Immunogen:	von Hippel-Lindau binding protein 1
Synonyms:	Prefoldin subunit 3 HIBBJ46 von Hippel-Lindau-binding protein 1 (VBP-1, VHL-binding protein 1) VBP1 PFDN3
Observed MW:	23-25 kDa
Uniprot ID :	P61758

Application

Reactivity:	Human, Mouse, Rat
Tested Application:	ELISA, IHC, IF, WB
Recommended dilution:	WB: 1:500-1:5000; IHC: 1:20-1:200;IF: 1:20-1:200
Image:	

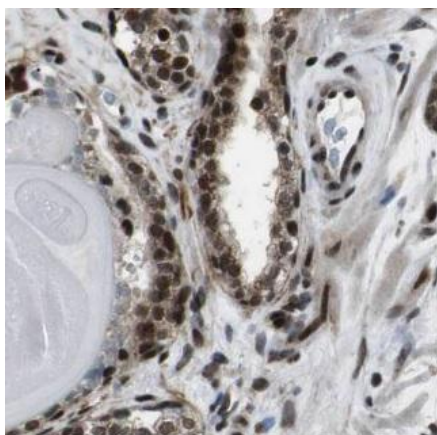

Immunohistochemistry of paraffin-embedded rat testis using FNab09379(VBP1 antibody) at dilution of 1:100

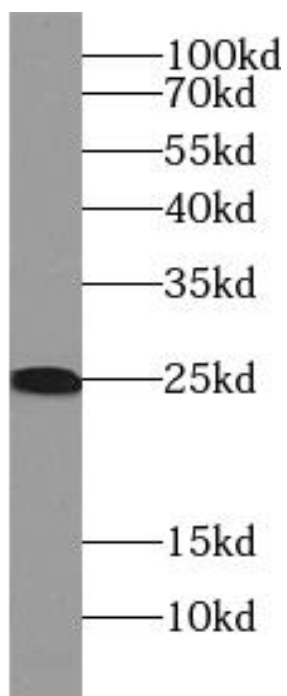

rat brain tissue were subjected to SDS PAGE followed by western blot with FNab09379(VBP1 Antibody) at dilution of 1:1000