

## ZNF35 antibody

### Product Information

Catalog No.:	FNab09689
Size:	100µg
Form:	liquid
Purification:	Immunogen affinity purified
Purity:	≥95% as determined by SDS-PAGE
Host:	Rabbit
Clonality:	polyclonal
Clone ID:	None
IsoType:	IgG
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3, -20°C for 12 months(Avoid repeated freeze / thaw cycles.)

### Background


May be involved in transcriptional regulation. Involved in cell differentiation and/or proliferation.

### Immunogen information

Immunogen:	zinc finger protein 35
Synonyms:	Zinc finger protein 35 Zinc finger protein HF.10 ZNF35
Observed MW:	59 kDa
Uniprot ID :	P13682

### Application

Reactivity:	Human, Mouse, Rat
Tested Application:	ELISA, WB
Recommended dilution: WB:	1:200-1:2000
Image:	

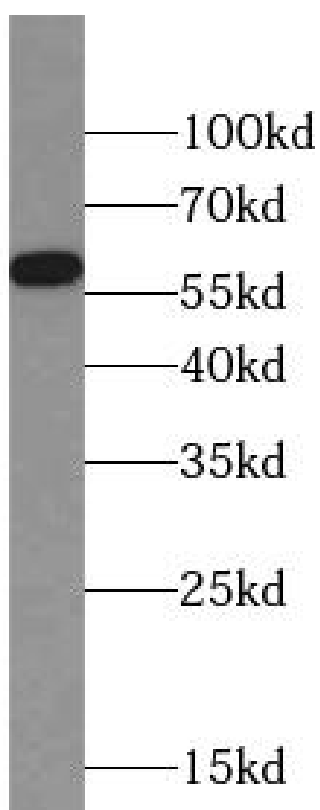

HeLa cells were subjected to SDS PAGE followed by western blot with FNab09689(ZNF35 antibody) at dilution of 1:300