

## PGAM1 antibody

### Product Information

Catalog No.:	FNab09874
Size:	100µg
Form:	liquid
Purification:	Immunogen affinity purified
Purity:	≥95% as determined by SDS-PAGE
Host:	Rabbit
Clonality:	polyclonal
Clone ID:	None
IsoType:	IgG
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3, -20°C for 12 months(Avoid repeated freeze / thaw cycles.)

### Background


Interconversion of 3-and 2-phosphoglycerate with 2,3-bisphosphoglycerate as the primer of the reaction. Can also catalyze the reaction of EC 5.4.2.4(synthase) and EC 3.1.3.13(phosphatase), but with a reduced activity.

### Immunogen information

Immunogen:	phosphoglycerate mutase 1(brain)
Synonyms:	Phosphoglycerate mutase 1 BPG-dependent PGAM 1 Phosphoglycerate mutase isozyme B (PGAM-B) PGAM1 PGAMA
Observed MW:	29 kDa
Uniprot ID :	P18669

### Application

Reactivity:	Human, Mouse, Rat
Tested Application:	ELISA, WB
Recommended dilution:WB:	1:200-1:2000
Image:	

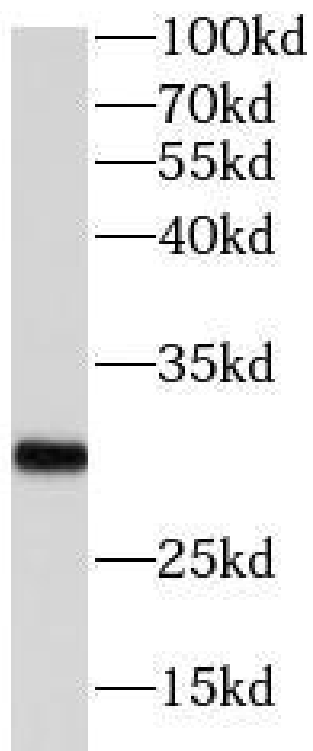

LO2 cells were subjected to SDS PAGE followed by western blot with FNab09874(PGAM1 antibody) at dilution of 1:1000