

anti- TAB3 antibody

Product Information

Catalog No.: FNab10439

Size: 100μg Form: liquid

Purification: Immunogen affinity purified

Purity: ≥95% as determined by SDS-PAGE

Host: Rabbit

Clonality: polyclonal

Clone ID: None IsoType: IgG

Storage: PBS with 0.02% sodium azide and 50% glycerol pH 7.3, -20°C for 12

months (Avoid repeated freeze / thaw cycles.)

Background

Adapter linking MAP3K7/TAK1 and TRAF6 or TRAF2. Mediator of MAP3K7 activation, respectively in the IL1 and TNF signaling pathways. Plays a role in activation of NF-kappa-B and AP1 transcription factor. Isoform 2 may be an oncogenic factor.

Immunogen information

Immunogen: TGF-beta-activated kinase 1 and MAP3K7-binding protein 3

Synonyms: Mitogen-activated protein kinase kinase kinase 7-interacting protein

3,NF-kappa-B-activating protein 1, TAK1-binding protein 3,

MAP3K7IP3

Observed MW: 75-100 kDa Uniprot ID: Q8N5C8

Application

Reactivity: Human, Mouse, Rat Tested Application: ELISA, WB, IHC, IF

Recommended dilution: WB: 1:500-1:1000; IHC: 1:20-1:200; IF: 1:20-1:100

Wuhan Fine Biotech Co., Ltd.


B9 Bld, High-Tech Medical Devices Park, No. 818 Gaoxin Ave. East Lake High-Tech Development Zone. Wuhan, Hubei, China (430206)

Tel:(0086)027-87384275 Fax: (0086)027-87800889 <u>www.fn-test.com</u>


1



Image:

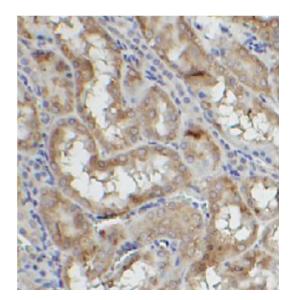

Immunohistochemistry of paraffin-embedded rat kidney tissue slide using FNab10439(TAB3 Antibody) at dilution of 1:50

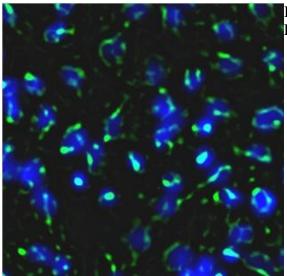


Immunofluorescent analysis of A431 cells using FNab10439 (TAB3 antibody) at dilution of 1:50.

Wuhan Fine Biotech Co., Ltd.

B9 Bld, High-Tech Medical Devices Park, No. 818 Gaoxin Ave. East Lake High-Tech Development Zone. Wuhan, Hubei, China (430206)

Tel:(0086)027-87384275 Fax: (0086)027-87800889 <u>www.fn-test.com</u>



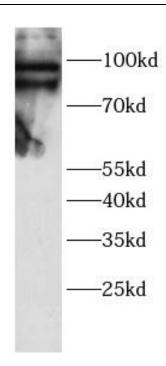

293T cells were subjected to SDS PAGE followed by western blot with FNab10439(TAB3 Antibody) at dilution of 1:1000