

---

## ATM antibody

### Product Information


Catalog No.:	FNab10667
Size:	100µg
Form:	liquid
Purification:	Immunogen affinity purified
Purity:	≥95% as determined by SDS-PAGE
Host:	Rabbit
Clonality:	polyclonal
IsoType:	IgG
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3, -20°C for 12 months (Avoid repeated freeze / thaw cycles.)

### Immunogen information

Immunogen:	ataxia telangiectasia mutated
Synonyms:	A T mutated, AT1, ATA, ataxia telangiectasia mutated, ATC, ATD, ATDC, ATE, ATM, Serine protein kinase ATM, TEL1, TELO1
Observed MW:	310-350 kDa
Uniprot ID :	Q13315

### Application

Reactivity:	Human
Tested Application:	ELISA, WB
Recommended dilution:WB:	1:500 - 1:1000
Image:	

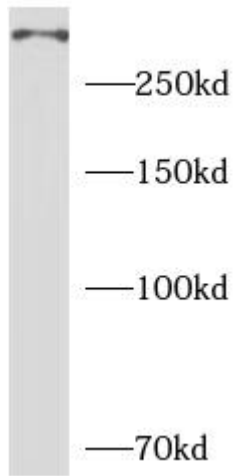

NCI-H460 were subjected to SDS PAGE followed by western blot with FNab10667(ATM antibody) at dilution of 1:500