

## PE Anti-Mouse CD127(A7R34)

### Product Information

Catalog No.:	PE-30013
Conjugation:	PE
Size:	50T/100T
Form:	liquid
Host:	Rat
Clonality:	monoclonal
Clone ID:	A7R34
IsoType:	IgG2a, kappa
Storage:	PBS with 0.1% sodium azide, 1%BSA, pH 7.3, 2-8°C for 12 months (Avoid repeated freeze / thaw cycles.)

### Immunogen information


Synonyms:	Interleukin-7 receptor subunit alpha IL7r IL-7 receptor subunit alpha IL-7R subunit alpha IL-7R-alpha IL-7RA CD127
Uniprot ID :	P16872

### Application

Reactivity:	Mouse
Tested Application:	FC

Recommended dilution: Volume per test: 5µL. Each lot of this antibody is quality control tested by flow cytometric analysis. The amount of the reagent is suggested to be used 5 µL of antibody per test (million cells in 100 µL staining volume or per 100 µL of whole blood). Please check your vial before the experiment. Since applications vary, the appropriate dilutions must be determined for individual use.

Image:

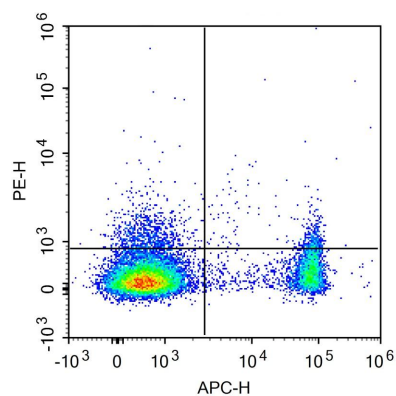

BALB/c murine splenocytes are stained with PE Anti-Mouse CD127 Antibody and APC Anti-Mouse CD4 Antibody.