


Recombinant Human CRYAB

Catalogue No.:	P1151
Species:	Human
Uniprot ID:	P02511
Expression Region:	1-175
Host:	E.Coli
Tags:	N-terminal His Tag
Molecular Weight:	22 kDa under reducing conditions
Purity:	Greater than 90% as determined by SDS-PAGE.
Formulation:	Lyophilized from a 0.2 µm filtered solution in 10 mM Hepes, 150 mM NaCl with 5% trehalose, pH7.4.
Reconstitution:	Centrifuge the vial before opening, reconstitute in sterile distilled water to a concentration of 0.1-1 mg/ml by gently pipetting 2-3 times, don't vortex.
Storage:	The lyophilized protein is stable at -20 °C for up to 1 year. After reconstitution, the protein solution is stable at -20 to -80 °C for 3 months or 1 week at 2 to 8 °C under sterile conditions. For extended storage, it is recommended to further dilute in working aliquots, avoid repeated freeze/thaw cycle.

Synonyms: P2RY11 antibody, P2Y purinoceptor 11 antibody, P2Y11 antibody, P2Y11 receptor antibody, P2Y11_HUMAN antibody, Purinergic receptor P2Y G protein coupled 11 antibody, Purinergic receptor P2Y11 antibody

SDS-PAGE:

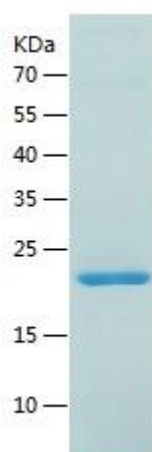

Safety Note: This product is intended for research and manufacturing uses only. It is not a diagnostic device. Product degradation will result from multiple freeze/thaw cycles. It is suggested that the antigen be stored in use size aliquots and thawed just prior to use. This material has been inactivated, however as with all biological materials, it should be handled as potentially infectious. The user assumes all responsibility for care, custody and control of the material, including its disposal, in accordance with all regulations.