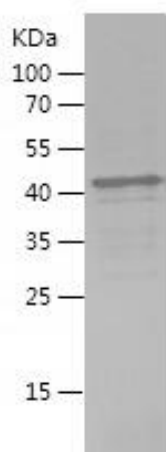


Recombinant Human WRN

Catalogue No.:	P1968
Species:	Human
Uniprot ID:	Q14191
Expression Region:	1224-1432
Host:	E.Coli
Tags:	N-terminal His-IF2DI Tag
Molecular Weight:	43.8 kDa under reducing conditions
Purity:	Greater than 90% as determined by SDS-PAGE.
Formulation:	Lyophilized from a 0.2 µm filtered solution in 10 mM Hepes, 150 mM NaCl with 5% trehalose, pH7.4.
Reconstitution:	Centrifuge the vial before opening, reconstitute in sterile distilled water to a concentration of 0.1-1 mg/ml by gently pipetting 2-3 times, don't vortex.
Storage:	The lyophilized protein is stable at -20 °C for up to 1 year. After reconstitution, the protein solution is stable at -20 to -80 °C for 3 months or 1 week at 2 to 8 °C under sterile conditions. For extended storage, it is recommended to further dilute in working aliquots, avoid repeated freeze/thaw cycle.

Synonyms:

DKFZp686C2056, DNA helicase, DNA helicase, RecQ like type 3, Exonuclease WRN, HGNC 12791, OTTHUMP00000225301, RecQ protein-like 2, RecQ-like type 3, RecQ3, RECQL2, RECQL3, Werner syndrome ATP-dependent helicase, Werner syndrome helicase, Werner syndrome protein, Werner syndrome, RecQ helicase like, WRN, WRN_HUMAN

SDS-PAGE:**Safety Note:**

This product is intended for research and manufacturing uses only. It is not a diagnostic device. Product degradation will result from multiple freeze/thaw cycles. It is suggested that the antigen be stored in use size aliquots and thawed just prior to use. This material has been inactivated, however as with all biological materials, it should be handled as potentially infectious. The user assumes all responsibility for care, custody and control of the material, including its disposal, in accordance with all regulations.