


Recombinant Human DYNLL1

Catalogue No.:	P5535
Species:	Human
Uniprot ID:	P63167
Expression Region:	1-89
Host:	E.Coli
Tags:	N-terminal His Tag
Molecular Weight:	9.7 kDa under reducing conditions
Purity:	Greater than 90% as determined by SDS-PAGE.
Formulation:	Lyophilized from a 0.2 µm filtered solution in 10 mM Hepes, 150 mM NaCl with 5% trehalose, pH7.4.
Reconstitution:	Centrifuge the vial before opening, reconstitute in sterile distilled water to a concentration of 0.1-1 mg/ml by gently pipetting 2-3 times, don't vortex.
Storage:	The lyophilized protein is stable at -20 °C for up to 1 year. After reconstitution, the protein solution is stable at -20 to -80 °C for 3 months or 1 week at 2 to 8 °C under sterile conditions. For extended storage, it is recommended to further dilute in working aliquots, avoid repeated freeze/thaw cycle.

Synonyms: 8 kDa dynein light chain, 8kDLC, Cytoplasmic dynein light polypeptide, DLC1, DLC8, DNCL1, DNCLC1, DYL1_HUMAN, Dynein , cytoplasmic, light chain 1, Dynein light chain 1 cytoplasmic, Dynein light chain 1, cytoplasmic, Dynein light chain LC8 type 1

SDS-PAGE:

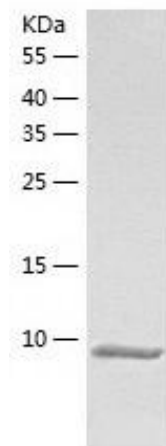

Safety Note: This product is intended for research and manufacturing uses only. It is not a diagnostic device. Product degradation will result from multiple freeze/thaw cycles. It is suggested that the antigen be stored in use size aliquots and thawed just prior to use. This material has been inactivated, however as with all biological materials, it should be handled as potentially infectious. The user assumes all responsibility for care, custody and control of the material, including its disposal, in accordance with all regulations.