


Recombinant Human PPA2

Catalogue No.:	P6658
Species:	Human
Uniprot ID:	Q9H2U2
Expression Region:	33-334
Host:	E.Coli
Tags:	N-terminal His Tag or N-terminal His-IF2DI Tag, determined during production process
Molecular Weight:	33.1 kDa under reducing conditions
Purity:	Greater than 90% as determined by SDS-PAGE.
Formulation:	Lyophilized from a 0.2 μ m filtered solution in 10 mM HEPES, 150 mM NaCl with 5% trehalose, pH7.4.
Reconstitution:	Centrifuge the vial before opening, reconstitute in sterile distilled water to a concentration of 0.1-1 mg/ml by gently pipetting 2-3 times, don't vortex.
Storage:	The lyophilized protein is stable at -20 °C for up to 1 year. After reconstitution, the protein solution is stable at -20 to -80 °C for 3 months or 1 week at 2 to 8 °C under sterile conditions. For extended storage, it is recommended to further dilute in working aliquots, avoid repeated freeze/thaw cycle.

Synonyms: HSPC124, Inorganic pyrophosphatase 2, Inorganic pyrophosphatase 2, mitochondrial, IPYR2_HUMAN, mitochondrial, Ppa2, PPase 2, PPase2, Pyrophosphatase (inorganic) 2, Pyrophosphatase SID6 306, Pyrophosphatase SID6-306, Pyrophosphate phospho-hydrolase 2, SID6 306

SDS-PAGE:

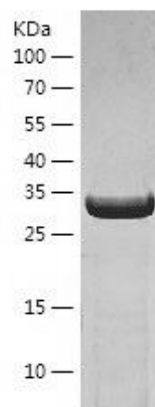

Safety Note: This product is intended for research and manufacturing uses only. It is not a diagnostic device. Product degradation will result from multiple freeze/thaw cycles. It is suggested that the antigen be stored in use size aliquots and thawed just prior to use. This material has been inactivated, however as with all biological materials, it should be handled as potentially infectious. The user assumes all responsibility for care, custody and control of the material, including its disposal, in accordance with all regulations.