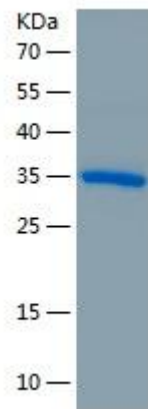


Recombinant Mouse CCL25

Catalogue No.:	P9665
Species:	Mouse
Uniprot ID:	O35903
Expression Region:	24-144
Host:	E.Coli
Tags:	N-terminal His-IF2DI Tag
Molecular Weight:	34.8 kDa under reducing conditions
Purity:	Greater than 95% as determined by SDS-PAGE
Formulation:	Lyophilized from a 0.2 μ m filtered solution in 10 mM Hepes, 500 mM NaCl with 5% trehalose, pH7.4.
Reconstitution:	Centrifuge the vial prior to opening, reconstitute in sterile distilled water or buffer of interest to a concentration of 0.1-1mg/ml by gently pipetting 2-3 times, don't vortex.
Storage:	The lyophilized protein is stable at -20°C for up to 1 year, the protein solution can be stored at 2-8°C for up to 1 week. For extended storage, it is recommended to further dilute in a buffer containing 0.1% BSA (make sure BSA will NOT interfere with your experiment) and store in working aliquots at -20°C to -80°C for 6 months. Avoid repeated freeze/thaw cycle.

Synonyms:

A130072A22Rik, C-C motif chemokine 25, CCL 25, CCL25, CCL25_HUMAN, Chemokine (C C motif) ligand 25, Chemokine (CC motif) ligand 25, Chemokine C C motif ligand 25, Chemokine CC motif ligand 25, Chemokine TECK, Ck beta 15, Ckb 15, Ckb15, MGC125074, MGC150327, SCY A25, SCYA 25, SCYA25, Small inducible cytokine A25

Image:**Safety note:**

This product is intended for research and manufacturing uses only. It is not a diagnostic device. Product degradation will result from multiple freeze/thaw cycles. It is suggested that the antigen be stored in use size aliquots and thawed just prior to use. This material has been inactivated, however as with all biological materials, it should be handled as potentially infectious. The user assumes all responsibility for care, custody and control of the material, including its disposal, in accordance with all regulations.